

nk medical products, inc.

Pads and Accessories



Choose from three different styles of O.R. Table Pads and Armboard Pads:

Standard

High-density polyurethane foam with a conductive vinyl sheeting top and bottom cover. This pad is recommended for shorter procedures.

Deluxe

High-density polyurethane foam with a soft, stretchable, urethane-coated, conductive top cover and a conductive vinyl sheeting bottom cover. This pad is recommended as an upgrade over the Standard Pad for short to medium length procedures.

Softcare

This pad incorporates a layer of pressure relief viscoelastic foam (a.k.a. Memory Foam) to aid in the prevention of pressure ulcers during longer procedures. The balance of the foam is high-density polyurethane foam. The top cover material is the same soft, stretchable, urethane-coated, conductive material that is used on the Deluxe Pad, with a conductive vinyl sheeting bottom cover.

- All cover materials are LATEX FREE, radiolucent, waterproof, stain resistant, anti-microbial, anti-bacterial, washable and conform to California Technical Bulletin 117 (CAL117).
- All foam fillers are LATEX FREE and conform to CAL 117.
- All completed pads conform to CAL 129.
- Boston Fire Code and/or CAL 603 compliance are available.




NK MEDICAL PRODUCTS, INC.
10123 MAIN STREET • CLARENCE, NEW YORK 14031
Toll Free: 1-800-274-2742 • (716) 759-7200 • Fax: (716) 759-0700
www.nkmedicalproducts.com ~ info@nkmedicalproducts.com




Common O.R. Pad Sets


AMSCO 1080/2080

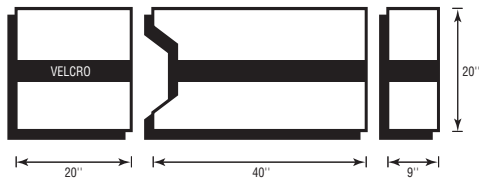


SIEMENS/MAQUET 1130

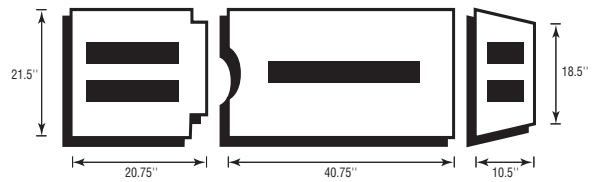


AMSCO 1080/2080 w/Foot Extension

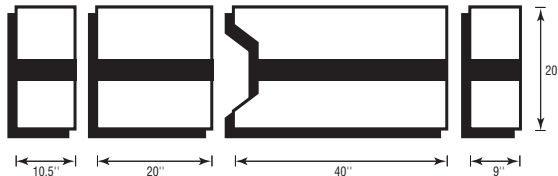


MIDMARK SOLID TABLE

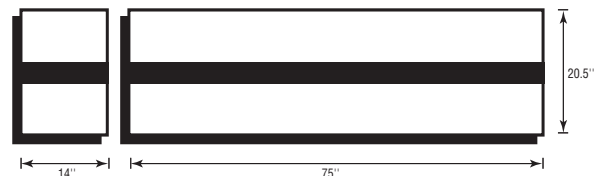


AMSCO 3080/3085

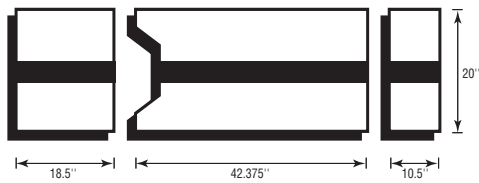


MIDMARK 7100

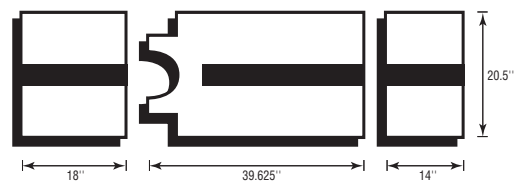


AMSCO URO-ENDO 1100

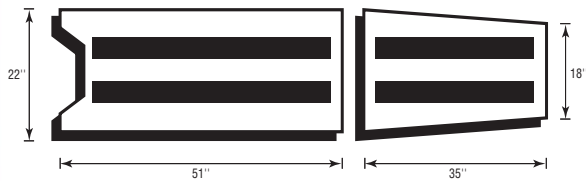


SKYTRON 3600

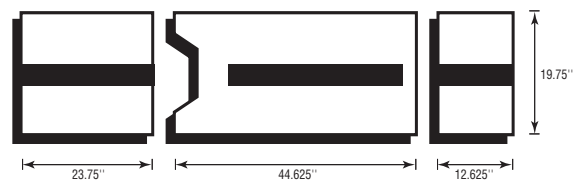


SHAMPAINE 1400/1500/1900

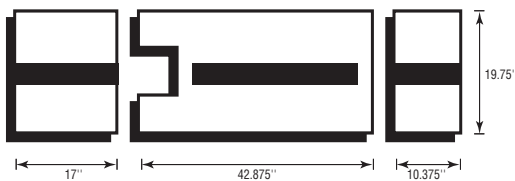

SKYTRON 3100

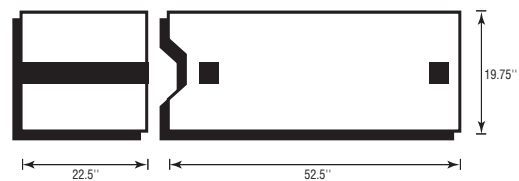


SHAMPAINE 4900B/5100B

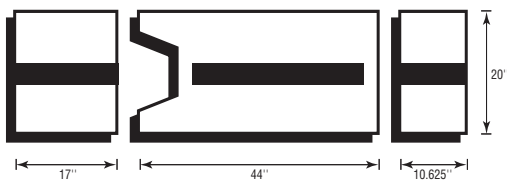


SKYTRON 3500



SHAMPAINE RAD-OP/1501

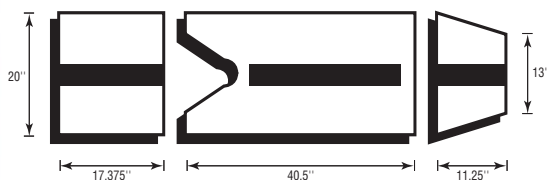


SKYTRON 5001/6001/6500/1201

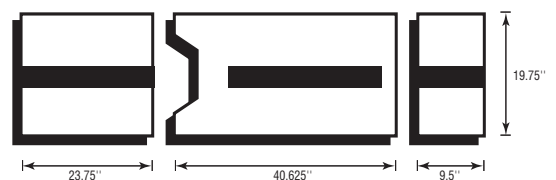


Common Stretcher/Gurney Pads


Stryker - Selection

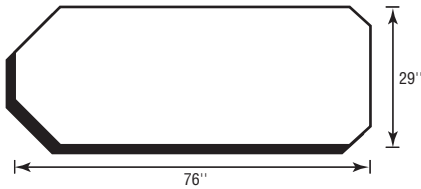


Hill-Rom - Selection

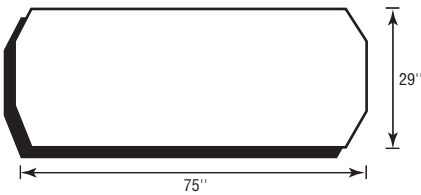


Hill-Rom Trans Star

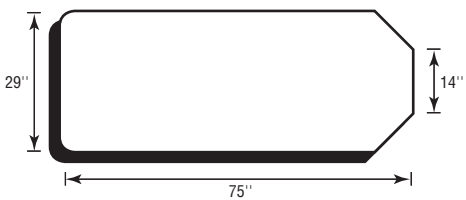

Hausted, Midmark and Others

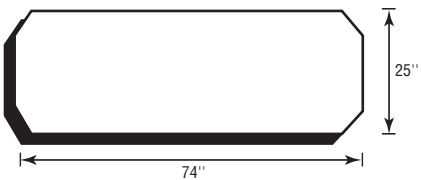

Universal - Fits All

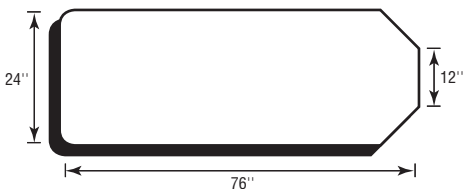

Choose from three different styles of Gurney / Stretcher Pads:

Standard

High-density polyurethane foam 3-4" thick with a durable, heavy-duty naugahyde top and bottom cover. Color variations and selected Department Names Decals are available on Standard Gurney Pads to help you keep track of "wandering" stretchers.

Deluxe

High-density polyurethane foam 3-4" thick with a soft, stretchable, urethane-coated, conductive top cover and a conductive vinyl sheeting bottom cover. This pad is recommended as an upgrade over the Standard Gurney Pad. Available in black only. Selected Department Name Decals are available.

Softcare

This 3-4" thick pad incorporates a layer of pressure relief viscoelastic foam (a.k.a. Memory Foam) to aid in the prevention of pressure ulcers. The balance of the foam is high-density polyurethane foam. The top cover material is the same soft, stretchable, urethane-coated, conductive material that is used on the Deluxe Pad, with a conductive vinyl sheeting bottom cover. Available in black only. Selected Department Name Decals are available.

- All cover materials are radiolucent, waterproof, stain resistant, anti-microbial, anti-bacterial, washable and conform to California Technical Bulletin 117 (CAL117).
- All foam fillers conform to CAL 117.
- Boston Fire Code and/or CAL 603 compliance are available.



Birthing Bed Mattresses



All Birthing Bed Mattresses are made with high-density polyurethane foam. Available in Standard Beige nylon material, or Comfort Blue light-stretch material.

Please specify Manufacturer's Name and Model Number when ordering. Many older birthing bed pads used snaps to adhere the pads to the table. Wherever possible, we will replace these snaps with Velcro to increase the durability of your pads.

- All cover materials are LATEX FREE, radiolucent, waterproof, stain resistant, anti-microbial, anti-bacterial, washable and conform to California Technical Bulletin 117 (CAL117).
- All foam fillers are LATEX FREE and conform to CAL 117.
- All completed pads conform to CAL 129.
- Boston Fire Code and/or CAL 603 compliance are available.

Restraint Straps

Conductive Patient Restraint Strap is 96" Long (enough length to wrap around the patient AND the table for tables without universal bars). Fits all surgical tables. Strong cam-buckles




lock strap tightly in place. Easy to clean, just wipe with soap and water or disinfectant. Universal Clips (shown in photo) are available for an extra charge – they snap onto the universal bar and STAY IN PLACE on their own. Extended length restraint straps are also available.

Arm Boards Pads


Choose from Standard, Deluxe or Softcare style arm board pads (see front cover for details). All Makes and Models are available.



AMSCO/Steris Armboard Pad

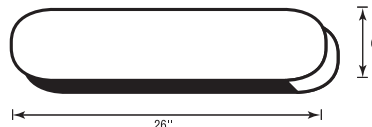

Skytron Armboard

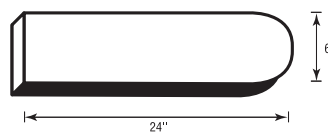

MDT Champagne and OSI Armboard Pad

Boston Fire Code Available

Tested and certified by an independent testing laboratory as meeting the Boston Fire Department IX-11 Mattress Fire Test. Pad contains a special internal barrier lining that shields the foam from fire.



NK MEDICAL PRODUCTS, INC.

10123 MAIN STREET • CLARENCE, NEW YORK 14031

Toll Free: 1-800-274-2742 • (716) 759-7200 • Fax: (716) 759-0700

www.nkmedicalproducts.com ~ info@nkmedicalproducts.com